

# Vela Tango 700 Configuration & Accessory Options

## Specifications

Frame size: 55 x 55cm

Chair's weight: 27kg (700), 34,5kg (700E)

Recommended user weight: max. 130kg (700), 145kg (700E)

Maximum weight load: 160kg

## Height Adjustment

The user's strength and energy levels during the day is a good indication to define how much help and support the user needs from the chair. You can choose between a manual (gas) lift, or electric lift.



**Electric: Vela Tango 700EL**



**Manual / Gas lift: Vela Tango 700**

### Electric height adjustment

Stepless, electric height adjustment of the seat provides gentle and safe regulation of the seat height.

The electric lift helps with standing up - with a touch of a button, the seat gently raises, and the user can safely stand up with a minimum effort. Or reach higher up cupboards, or reach lower and unload a dishwasher.

The button is placed under the armrest and can be fitted on the left or right side.

## Build height

When choosing a seat height, it is important to consider the user's size, level of function, needs and where the chair will be used. Note: if your height is below 165cm, we recommend **Low** seat height.

### Electric: Vela Tango 700EL

(user weight max. 145kg)

	Sitting height with Active seat	Sitting height with Contour seat
Low	41-61cm	43-63cm
Medium	47.5 - 77.5cm	49.5 - 79.5cm
High	55 - 85cm	57-87cm

### Manual / Gas lift: Vela Tango 700

(user weight max. 130kg)

	Sitting height with Active seat	Sitting height with Contour seat
Low	38.5 - 48.5cm	40.5 - 50.5cm
Medium	46.5 - 61.5cm	48.5 - 63.5cm
High	55 - 75cm	57-77cm

## Seats

When choosing a seat type, it is important to know the user's size and support needs when sitting down.

Adjustable seat angle: 10° forwards, 15° backwards.

Seat rotation: 90° to the right/90° to the left



### Active Seat

Offers great comfort and light support. It comes in various sizes to fit the user's body and support needs.

#### Available in:

Small: W40 x D40cm  
Medium: W44 x D44cm  
Large: W47 x D47cm  
XLarge: W49 x D49cm  
Deep: W40 x D49cm



### Contour Seat

Offers excellent comfort and a high degree of support. The seat is shaped with side support that helps stabilise your sitting position.

#### Available in:

Small: W40 x D40cm  
Medium: W44 x D44cm  
Large: W47 x D47cm  
XLarge: W49 x D49cm



### Seatplate only

This allows for the user to use a custom cushion like their wheelchair cushion or a pressure care cushion on the chair.

#### Available in:

Small: W40 x D48cm  
Medium: W45 x D48cm  
Large: W50 x D48cm

## Backrests

When choosing a backrest, be aware of the user's need for support versus free movement. A backrest shall provide support and stability at the lumbar and upper part of the pelvis so that the user receives support in the sitting position.

Adjustable backrest angle: 4° forwards, 22° backward

Height adjustment of backrest: 10 cm



### Active Backrest

Offers great comfort and light support.

#### Available in:

Small: W38 x H22cm  
Medium: W38 x H30cm  
Large: W40 x H40cm  
X Large: W40 x H48cm




### Contour Backrest

Offers excellent comfort and a high degree of support. The backrest has contours with side support that helps stabilise your sitting position.

#### Available in:

Medium: W36 x H38cm  
Large: W39 x H43cm  
X Large: W41 x H48cm

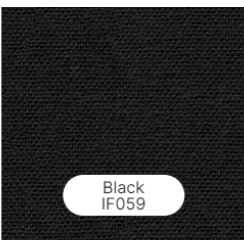

## Fabric

All Vela chairs from Morton & Perry are upholstered with 'Main Line Plus' fabric in the colour black. Main Line Plus is a popular textile of high quality from Camira made from 67% Wool, 20% Flame-Retardant Viscose and 13% Viscose. With a wear resistance of >50,000 Martindale, it is known for having a great resistance to wear, good comfort and high cleanability. Other fabrics are available on indent orders only. The upholstery is cleaned by vacuuming, professional cleaning or wiping with a clean, wrung out cloth. A seat cover is available to order.

## Armrests

The **adjustable armrest** can be adjusted in height and width.

The **foldable armrest** has the same adjustment features as the 'adjustable armrest' but in addition is easily folded down, e.g. for sideways transfer, or help with getting dressed.

Cushion material: Soft cushion made from PU-padding

Dimensions: width x length 9 x 25 cm

Armrest height: adjustable 14-25 cm

Distance between armrests: adjustable 38,5-52,5cm



### Armrest adjustable

Adjust armrest height & width.



### Armrest foldable

Adjust armrest height, width and drop back armrest for easy transfers.

---

## Brake

A reliable brake is an essential feature on the VELA chair, as the brake creates the necessary stability and safety that makes activity, learning and independence possible.



### Brakerod standard 39cm

This is the chair's standard brake. We recommend it if the person has no special needs, and the chair is at a standard height.



### Brakerod Adjustable 44-62cm

The extension handle extends the brake lever, making it easier to reach, compared to a standard brake. We recommend the extension handle if the person is mainly sitting up high. And extension handle also makes it easier to brake the chair, as you have to use less force than with a standard brake handle.



### Brakerod Adjustable 44-62cm with T-lever

The T-handle extends the brake lever, and at the same time, makes it wider, so it is easier to reach. We recommend this handle for people who have reduced mobility in the shoulders, arms or upper back. The T-handle makes it easier to brake the chair, as you have to use less force than with a standard handle.

## Footrests

A footrest is often essential for the user to have an optimal sitting position. The feet must be supported when the user is sitting in the chair, so if the user cannot reach the floor while sitting, a footrest must be used.



### Footrest fold-up with legs

We recommend a fold-up footrest if the person sometimes sits at a height where it is impossible to rest their feet on the floor.



### Footrest swing-away set

We recommend single footrests for one or both feet if you need light support and often switch between high and low sitting height. The footrests can be rotated with the foot and provide support and rotated away again when the user needs to move the chair forward.

---

## Anti-tip bar


We recommend an anti-tip bar on all chairs for users who have involuntary movements or uses the chair with a pedal exerciser. A 4-wheel brake is a benefit in those cases.



---

## Push bar adjustable

A push bar optimises the carers working position when there is a need to transport a person in the chair over short distances.

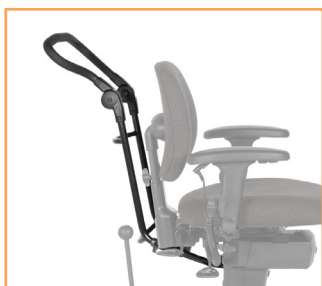

This push bar can be easily adjusted in angle, height, and depth to suit the person who is pushing the chair. The push bar is easy to remove. We recommend a push bar if a carer needs to transport a person, for example, to another room.

# Vela Tango 700 Prescription order form

Name: \_\_\_\_\_ Phone: \_\_\_\_\_ Email: \_\_\_\_\_

Address: \_\_\_\_\_

## Seat *\*Compulsory. Only tick one box please*

### Active Seat

- Small: W40 x D40cm
- Medium: W44 x D44cm
- Large: W47 x D47cm
- X Large: W49 x D49cm
- Deep: W40 x D49cm

### Contour Seat

- Small: W40 x D40cm
- Medium: W44 x D44cm
- Large: W47 x D47cm
- X Large: W49 x D49cm

### Seatplate only

- Small: W40 x D48cm
- Medium: W45 x D48cm
- Large: W50 x D48cm

## Backrest *\*Compulsory. Only tick one box please*

### Active Backrest

- Small: W38 x H22cm
- Medium: W38 x H30cm
- Large: W40 x H40cm
- X Large: W40 x H48cm

### Contour Backrest

- Medium: W36 x H38cm
- Large: W39 x H43cm
- X Large: W41 x H48cm

## Lift & Build Height *\*Compulsory. Choose either Electric or Gas lift. Then select: Low, Medium or High.*

- Electric: Vela Tango 700EL**  
(user weight max. 145kg)

	Sitting height with Active seat	Sitting height with Contour seat
<input type="checkbox"/> Low	41-61cm	43-63cm
<input type="checkbox"/> Medium	47.5 - 77.5cm	49.5 - 79.5cm
<input type="checkbox"/> High	55 - 85cm	57-87cm

- Manual / Gas lift: Vela Tango 700**  
(user weight max. 130kg)

	Sitting height with Active seat	Sitting height with Contour seat
<input type="checkbox"/> Low	38.5 - 48.5cm	40.5 - 50.5cm
<input type="checkbox"/> Medium	46.5 - 61.5cm	48.5 - 63.5cm
<input type="checkbox"/> High	55 - 75cm	57-77cm

## Armrests *\*Compulsory. Only tick one box please*

- Standard  Foldable

## Electric Elevate button (under armrest)

- Right hand side **or**  Left hand side

## Brake *\*Compulsory. Only tick one box please*

- Standard Brakerod 39cm  Brakerod adjustable 44-62cm  Brakerod adjustable 44-62cm with T-Lever

## Footrest *Optional. Only tick one box please*

- Footrest fold-up with legs  Footrest swing-away set

## Accessories *Optional. Tick one or both*

- Anti-tip bar  Pushbar adjustable